

# PRODUKTDATENBLATT



<b>Produkt</b>	FFP2 UNICO3 CAP Atemschutzmaske (Einweg)
<b>Artikelnr.</b>	cap-med-013
<b>Verpackungseinheit (VPE)</b>	10/10 Stk. Box / 1.350 Stk. im Karton (Verpackung in Englisch)
<b>Beschreibung</b>	FFP2 Einweg Atemschutzmaske ohne Filter, flach faltbar in der Mitte, mit elastischen Ohrschlaufen
<b>Farbe</b>	Weiß
<b>Zertifizierung</b>	DIN CE/EN 149:2001 + A1:2009 FFP NR D* CE Body Notified Identifikations Nr. 2163 *D = Sie haben den Dolomitenstaubtest bestanden und sind daher geeignet für die Verwendung in einer sehr staubigen Umgebung
<b>Schutzklasse</b>	FFP2, BFE 95% hohe Bakterienfilterung
<b>Anwendungsbereich</b>	Bakterien und Viren, luftgetragene Stäube, Allergene, Sporen und Hautschuppen

CAPULUS übernimmt keine Gewährleistung/ Haftung für die Aktualität, Richtigkeit und Vollständigkeit der vom Hersteller bereitgestellten Informationen und Daten. Das Produkt kann vom o.g. Bild abweichen.

## Kontakt / Bestellungen:

Tel: +49 40 307 772 15

Zertifizierung >> hier anklicken <<

E-Mail: info@capulus-medical.de


Anfragen: Bestellformular Geschäftskunden

Internet: www.capulus-medical.de

# FILTERING

## CASSTAR

### SMS Five-tier Structure

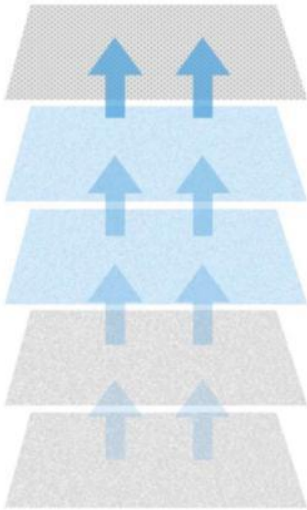

- **Non-woven surface :**  
Block larger particles;
- **High-efficiency filter layer :**  
Effective filtering of particles  $\geq 0.075$  microns
- **High-efficiency filter layer :**  
Effective filtering of particles  $> 0.075$  microns
- **Non-woven inner layer :**  
Skin friendly, comfortable and breathable
- **Non-woven inner layer :**  
Skin friendly, comfortable and breathable

## CASSTAR

### 3D Stereo Surround Design




Low expiratory resistance ,efficient filtration  
fluid resistant to splash and spatter of  
blood and other infectious material

## CASSTAR


### Five Layer Filtration

Blocking harmful substances layer by layer

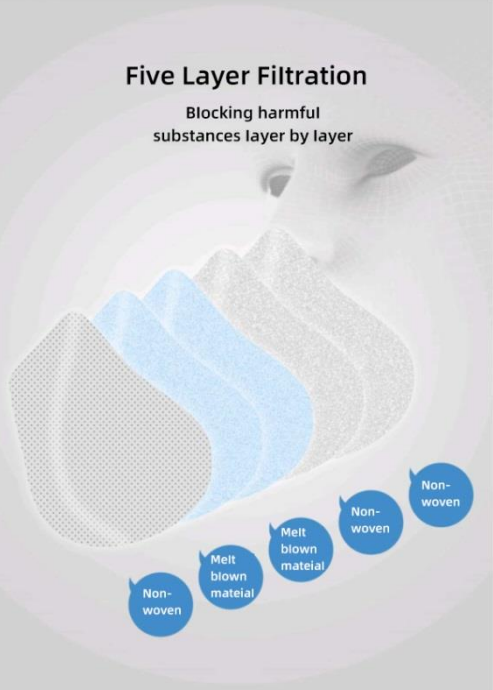

# ZERTIFIZIERUNG

UNIVERSAL



**UNIVERSAL**  
CERTIFICATION

**NB 2163**

## CERTIFICATE OF CONFORMANCE

**Certificate No: 2163-PPE-1162/01**

Respiratory protective devices, filtering half masks to protect against particles manufactured by  
  
 Third Floor, Building 4, No. 87, Ulusoyun, South Road, New City, Halki Street, Fongze

Continues to fulfil the requirements of

**EN 149:2001 + A1:2009 Respiratory Protective Devices -  
Filtering Half Masks to Protect Against Particles -  
Requirements, Testing, Marking**

Based on the evaluation of test reports and internal quality control audit reports according to EN 149+A1:2009 and Personal Protective Equipment Regulation (EU) 2016/425 Annex VII (Module C2). This certificate implies that the manufactured products show below are in conformance with the approved EU Type Examination model and meets the requirements of the regulation.


**Product Definition**

Model	Class	EU Type Examination Certificate		
		Serial No	Date	Issuing NB No
TY-95	FFP2 NR	2163-PPE-1162	27.07.2020	2163


Here by the manufacturer is allowed to use notified body number (2163) and can fix CE mark, as shown below, on the Category III product models given above, with:

- Issuing an appropriate EU Declaration of Conformity according to **Personal Protective Equipment Regulation (EU) 2016/425 Annex 9**.
- Taking all measures necessary so that the manufacturing process and its monitoring ensure the homogeneity of production and conformity of the manufactured PPE with the type described in the EU type examination certificate.

This certificate is issued on **27/07/2020** and will be valid for one year, until **26/07/2021** if the manufacturer makes no major change in the product designs and manufacturing processes affecting the product performance on the essential health and safety requirement.




2163



Suat KACMAZ  
UNIVERSAL CERTIFICATION  
Director

Verify the validity with the QR code



Necip Fazil Bulvarı Keyap Sitesi E2 Blok No:44/84 Yukarı Dudullu Ümraniye - İSTANBUL - TÜRKYE - T:+90 216 455 80 80

[UNIVERSALCERT.COM](http://UNIVERSALCERT.COM)

